

Vierra **LiFe**

Lithium powered scooters

**LiFetime Battery
Guarantee⁺**



RASCAL
LiFe **BRITISH
MADE**

The Lithium Ferro-Phosphate Battery
is **'built to last'** & the **SAFEST**
Lithium Battery in the Mobility sector

RASCAL 

Rascal Vierra **LiFe** Profile:

The **Rascal Vierra LiFe** is a quality engineered, lightweight, Aluminium boot scooter with State-Of-The-Art **Lithium Ferro-Phosphate Battery**. With a choice of 25Ah or 12.4Ah battery, this compact, class-leading Rascal includes a **Lifetime Battery Guarantee** for peace of mind. If you are looking for spacious legroom, all-round suspension with exceptional ride comfort, stability & optimum longevity, with this 4mph scooter offering upto 21.4 stone user weight, it is the obvious choice!

Greater flexibility with two easy access charging points, with **super light front section weighing just 10.8kg**, it is designed to dismantle with ease, this robust & highly manoeuvrable travel scooter with tight turning circle is great for both indoor & outdoor use.

Benefit from a range of upto 26.2 miles with the **Rascal Vierra LiFe**. Whatever journey you embark on, feel the level of comfort from the deluxe seat & enjoy the convenience of the smooth user-friendly controls, padded adjustable delta tiller & handy front basket which is quick to attach on or off for daily essentials.

LiFe moves you in a new direction...

The enhanced safety features include front & rear integrated **LED** high performance **lighting & Rapid 5A LiFePO4 charger** (available on 25Ah battery) providing upto 5 miles range per hour of charging.

Thoughtful practicality is packed into this modern scooter for both battery options. With the 12.4Ah battery, giving upto 13 miles range.



Discover your style with a choice of 5 vibrant colours, relax knowing that your stylish scooter will stand up to everyday life.




† LiFetime Battery Guarantee T&Cs apply

The **Rascal Vierra LiFe** with high quality State-Of-The-Art **Lithium Ferro-Phosphate** (LiFePO₄) Battery. Relax & enjoy the ride in the knowledge that the Lithium Ferro-Phosphate (LiFePO₄) Battery pack is the safest, **greenest** & is warranted for a design life of over 1500 cycles providing an increased efficiency & performance. Plus the LiFetime Battery Guarantee[†].

Main Features:

- British Made / Medical /Military grade **Lithium Ferro-Phosphate Battery**
- Lightweight, aluminium, high quality rugged frame
- All-round Suspension with exceptional ride comfort
- Front & Rear integrated LED high performance lighting
- 1150mm turning radius for outstanding manoeuvrability
- Tiller locking for easy transportation
- Fold down backrest with handy pouch & simple to remove adjustable armrests
- Front basket has good size & depth
- Two easy access charging points on tiller & battery box
- 8" Dual compound puncture proof tyres for impressive stability
- Folded Dimensions (LxWxH) 108cm x 47cm x 71cm / 42.6" x 18.6" x 27.9"
- Airline Compliant 12.4Ah battery weighs just 4.1kgs / 9 lbs

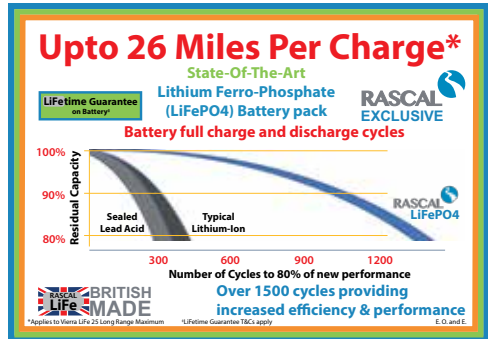

LiFetime Battery Guarantee

Safely transport your belongings & travel essentials with this compact accessory

OPTIONAL



Practical Features:



Air travel & CAA Lithium Battery Legislation
Airline Compliant 12.4Ah battery meets I.A.T.A. UN38.3 Passenger Aircraft Transportation requirements for overseas travel & can be taken on a plane. (Subject to Airline Regulation & Restrictions)

The **RASCAL LiFe** **BRITISH MADE** Lithium Ferro-Phosphate Battery is **'built to last'** & the **SAFEST** Lithium Battery in the Mobility sector

Specification:

	Vierra LiFe25Ah with Rapid 5A Charger	Vierra LiFe12.4Ah with 2A Charger
Type Class	B / 2	✓
Overall length	108.3cm / 42.6"	✓
Overall width	47cm / 18.6"	✓
Overall height	88cm / 34.7"	✓
Floor to seat height	39-43cm (min-max)	✓
Seat width	44cm / 17.3"	44cm / 17.3"
Seat depth	41cm / 16.1"	41cm / 16.1"
Folded Dimensions LxWxH	108cm / 47cm / 71cm (42.6" / 18.6" / 27.9")	✓
Maximum carrying capacity	136kgs / 21.4 stone	✓
Battery voltage and capacity	25.6V / 25Ah	25.6V / 6.2Ah x 2
Maximum speed	6.4km/h / 4mph	✓
Ground clearance	4.5cm / 1.8"	✓
Maximum kerb height	5cm / 2"	✓
Range*	up to 42.1km / 26.2 miles	up to 21km / 13 miles
The maximum safe slope	8°	✓
Turn-around width	1.25m / 4'1"	✓
Colour	Blue / Green / Orange / Pink / Yellow	
Single Heaviest component (front section)	10.8kgs / 23.8lbs	✓
Product weight (Excluding battery)	29.7kgs / 65.5 lbs	✓
(Including battery)	37.7kgs / 83.1 lbs	33.8kgs / 74.5 lbs

*Range on full charge and flat ground

See full Tech Page on Website & Owner's Manual for details

Explore with Confidence

Local Rascal Retailer:



† LiFetime Battery Guarantee T&Cs apply