

# Iconic

CosiChair

Tilt-in-Space Quad-Motor Riser Recliner



## Superior & Sophisticated

- Fully adjustable powered head & lumbar support
- True footrest extension
- Longest hip to heel seat
- Fitted USB & C charging ports



Featuring Zero Gravity seating position



# Iconic Quad-Motor Riser Recliner

Lasting relaxation, feel the difference



The Iconic is a Premium Riser Recliner designed by **CosiChair** which has a unique Tilt-in-Space (TIS) Quad-Motor function featuring the Zero Gravity position which evenly distributes the force & stress across the body, giving supported comfort from head to toe throughout movement offering greater relaxation.

Enjoy the adjustable powered head & lumbar positioning which offers personalised support for the head, neck & lower back. Benefit & take advantage of these features to watch TV, read a book or sit up to have a drink.

Designed using the finest fabrics & handmade to the highest standards, the Iconic offers first class flexibility with regards to **movement** & positioning, with Variable Angle Lift (VAL) adjusting the tilt angle on lift & the distinctive, automatic adjusting backrest opens up the angle for exit. This precision & smart technology allows a safe, controlled standing position with the user-friendly 9 button handset.



Blackberry



Natural



Smoke



Spring



Teal



Thunderstorm

## Features & Benefits


- Zero Gravity Seating Position
- Variable Angle Lift (VAL) to adjust the tilt angle on the lift function between 30°, 20° or 10° which reduces the seated angle when lifting & exiting the chair
- Very high leg lift giving full support behind the knee helping to relax the legs
- A 'Reset' button which has software designed to bring you back to a seated position in a safe & comfortable controlled manner from a recline or lift position
- Unique automatic adjusting backrest which helps to reduce the back angle to allow a more controlled safe standing position on exiting the chair
- Longest Hip to Heel seating with an additional 3" true footrest extension
- Contemporary Lateral or Waterfall Back Design
- Industry leading 7 Year Guarantee – T&C's apply\*



Available to complement your home in six modern, stylish & soft fabrics. Choose from the Lateral or Waterfall sumptuous back design or indulge yourself & have both...

Other designs, upholstery & colourways to suit all tastes are available on request.



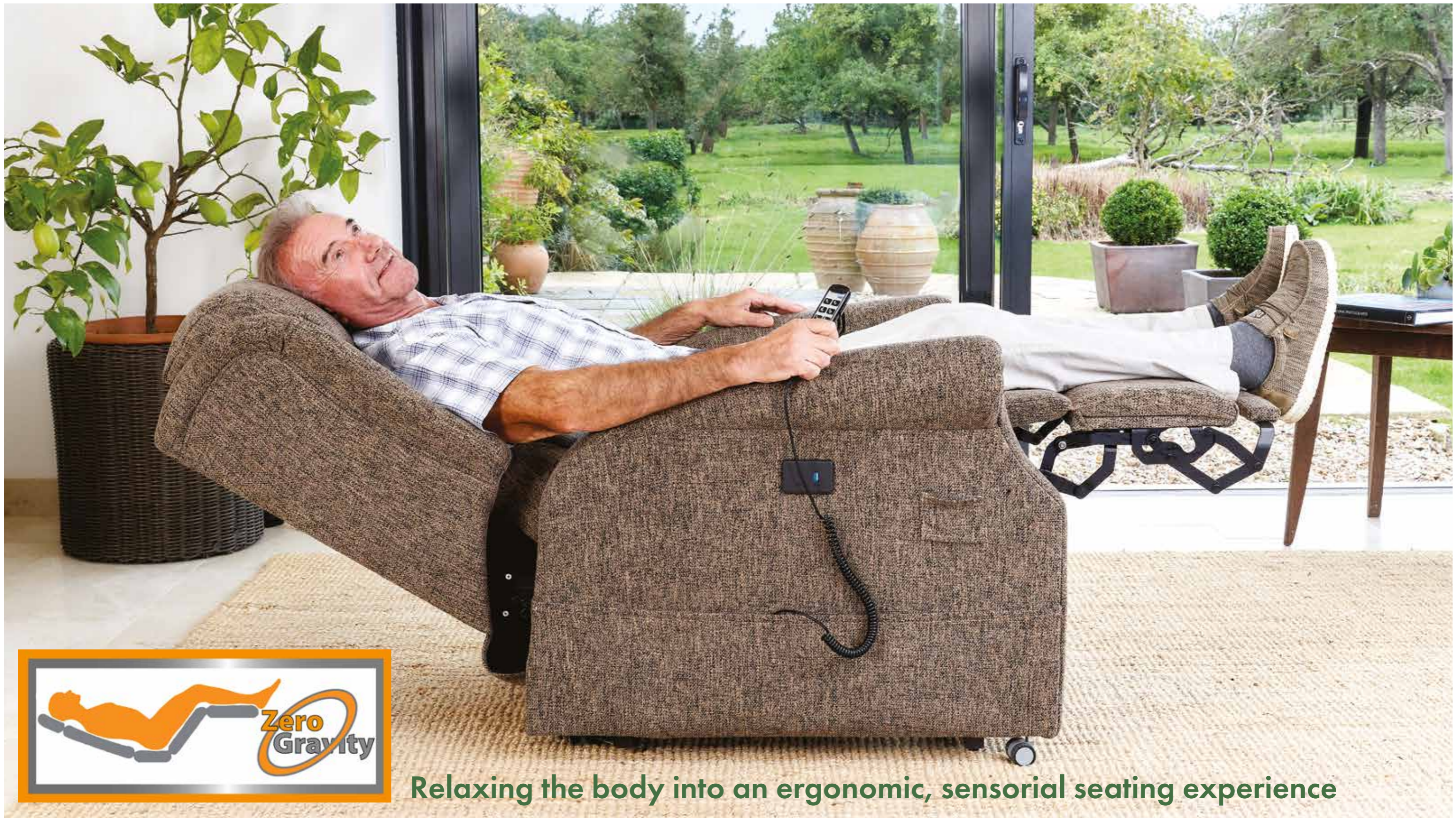

Relaxing the body into an ergonomic, sensorial seating experience



## Zero Gravity

The Zero gravity position puts your body into “zero stress zone” which offers even weight distribution across the body. Your muscles relax aiding breathing & as your legs are elevated this reduces stress on the heart. The position spreads your body weight out, reducing the stress & tension on the spine, hips, knees & lower back, benefitting your health & wellbeing.



Personalised support for the head, neck & lower back

### Truly improve your well-being...



A true extending footrest extension which supports the heel reducing the pressure placed upon the achilles



Fully adjustable powered headrest & neck support moves upto 14cm / 5½” for a perfect position

Fully adjustable powered lower lumbar support moves upto 5.5cm / 2 ¼ inches

\*7 Year Guarantee mechanism & electronics warranty. Handset & transformer fair wear & tear & a lifetime wooden frame warranty. We reserve the right to amend the duration of Component Guarantee Periods identified.

## Variable Angle Lift

10 degree

20 degree

30 degree



### Are you sitting comfortably?

The Iconic Riser Recliner has a Variable Angle Lift (VAL). There are 3 angles at which the chair can be set to rise at. If it isn't comfortable rising at the standard 30 degree angle, it can be set to a medium rise of 20 degrees or a “flat rise” which is just 10 degrees. The varying degrees of lift may help people that have postural problems, allowing supported & smoother transition from sitting to standing. There is the option to configure the VAL at time of purchase. Ask your Cosi Chair Specialist if you have specific requirements.

### The Revolutionary Iconic Quad-Motor comes with the following accessories as standard:

9v battery back up transformer

Side arm mounted USB & C ports which allows you to plug in & charge your electronic devices  
Twin pockets



Antimacassar head cover

Full length arm covers

Twin handset holder for left or right-handed control

Braked castors

**ELECTRIC  
MOBILITY**<sup>®</sup>



**Choose our furniture with confidence**

**Iconic Size Guide:**

	Small	Medium	Large
Seat height	47cm / 18½"	50.8cm / 20"	53.3cm / 21"
Seat width	48.3cm / 19"	50.8cm / 20"	53.3cm / 21"
Seat depth	45.7cm / 18"	50.8cm / 20"	55.9cm / 22"
Arm height	17.8cm / 7"	19.1cm / 7½"	19.1cm / 7½"
Back height	68.6cm / 27"	71.1cm / 28"	73.7cm / 29"
Hip to Heel – raised footrest	101.6cm / 40"	108cm / 42½"	114.3cm / 45"
Zero Gravity position	175.2cm / 69"	182.9cm / 72"	190.5cm / 75"
Footrest extension	7.6cm / 3"	7.6cm / 3"	7.6cm / 3"
Overall width	76.2cm / 30"	81.3cm / 32"	86.4cm / 34"
Overall height	91.4cm / 36"	94cm / 37"	96.5cm / 38"
Weight limit	25 stone / 158.8kg	25 stone / 158.8kg	25 stone / 158.8kg

Local Cosi Chair Specialist



Exclusively available from your local Cosi Chair Specialist

Due to a policy of continual improvements, Electric Mobility Euro Ltd reserve the right to change specification without prior notice. Images are for illustrative purposes only. Details are correct at time of going to print.

E. and O.E- v3 May 2024